

The precise FreeFold design delivers a smooth fold that requires 40% less effort than its closest competitor, ensuring a smooth folding experience with significantly less effort.

The cross-brace centers over the frame, making it as narrow as possible when folded. With a usable transport weight\* of just 16.1 lbs, getting it into your car is easier than ever.

Crafted with Sunrise Medical's proprietary technologies like specially heat-treated aluminum, ovalized cross-braces, and forged parts, ensuring an unparalleled riding experience.

The QS5 X features a strong box frame and ovalized 4-arm crossbrace, ensuring durability. Its sleek, modern design adds a stylish and sophisticated touch.


With two frame versions, multiple frame angles, and clinical seating options, the QS5 X ensures adaptability to diverse individual needs. This applies to the first-time rider or an experienced wheelchair rider, accommodating changes over time.


The optional no charge FOLD-N-LOCK feature conveniently keeps the chair securely folded with a simple "click". There's no need to hold the frame closed while lifting or storing.

Personalize your chair with a choice of frame graphics, painted handrims, painted cross-brace, accent upholstery trim, anodized components, and 25 paint colors.

Our online visualizer will help you to create the QS5 X of your dreams!



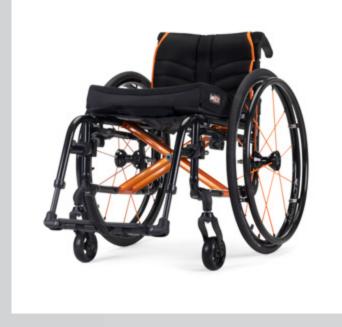

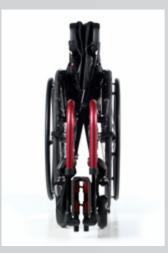



TECHNICAL DATA	
HCPCS CODE	K0005 / E1236
FRAME MATERIAL	7000 Aluminum
SEAT WIDTH	13" to 20"
SEAT DEPTH	13" to 20"
FRAME ANGLES	Swing Away: 60/70/80°, Fixed-Front: 75/85° and 1" inset
SEAT HEIGHT FRONT	14" to 22"
SEAT HEIGHT REAR	14" to 20"
BACK HEIGHT ADJUSTMENT	10" to 20"
BACK ANGLE ADJUSTMENT	From 80° to 105°
CENTER OF GRAVITY (COG)	1" to 3"
CHAIR WEIGHT *	Swing-Away: from 22.3 lbs., Fixed-Front: from 21.8 lbs.
USABLE TRANSPORT WEIGHT*	Swing-Away: from 16.5 lbs., Fixed-Front: from 16.1 lbs.
WEIGHT CAPACITY	300 lbs.
TRANSIT	Yes **

Refer to order form as specifications may vary based on configuration.

\* Chair weight is calculated based on a 16"x16" seat dimension with back upholstery, side guards & wheel locks. Transport weight is calculated with items that can be easily removed (i.e. Swing-Away footrests, armrests, rear wheels, and seating products).

Learn more



\*\* Conforms to ISO 7176-19 Transit Approved standards

