



ZIPPIE FAMILY

When you think of ZIPPIE, **imagine all the activities children can enjoy**. Each ZIPPIE product is designed to be a **simple-to-use** and **reliable support** for the lifestyles of children. ZIPPIE's superior designs increase **functional mobility** to keep children going, and our uncompromising quality ensures we **never let them down**.



ZIPPIE VOYAGE



QUICK FOLDING,
LIGHTWEIGHT,
STROLLER FRAME
WITH INCREDIBLY
VERSATILE SEATING
AND POSITIONING,
GIVING MORE OPTIONS
TO SUPPORT GROWING
BABIES AND THEIR
CAREGIVERS



- Voyage Base
- One-step quick-fold technology
 - Angle adjustable stroller handle
 - Airless wheels
 - Swivel front wheels that lock-out for long-distance strolling
 - Transit tie-downs



- Slide-N-Lock®
Quickly and comfortably attach and remove the seat from the stroller in one simple motion.

Many Functional
Accessories Available

Versatile Seating for Babies on-the-Grow!

The Voyage features a 35° tilt range and optional recline for superior positioning. The reversible seat can face forward so the child can explore their surroundings, or rearward for maximum bonding time. Multiple head, hip, thigh, and trunk support options create a personal fit that adjusts to accommodate rapid growth.



Quick Recline

Recline can be adjusted 85° to 100° whenever needed.



Convenient Tilt

The static tilt can be adjusted to 10°, 20° or 30° to best accommodate the child's needs.



The Most Compact Transport Dimensions in its class!

Transport dimensions (34" x 25.5" x 16.75") are small enough to fit within a standard compact car's trunk space, and even within some suitcases.

THE FIRST OF
ITS KIND: A TRULY
TRANSPORTABLE
PEDIATRIC TILT
WHEELCHAIR,
SPECIALLY DESIGNED
FOR FAMILIES
ON THE MOVE

ZIPPIE SPHYNX



ZIPPIE TS



TILT-IN-SPACE WHEELCHAIR IS AVAILABLE WITH A WIDE VARIETY OF OPTIONS, INCLUDING AN OPTIONAL FOLDING FRAME FOR EASY STORAGE.



45° of Tilt-in-Space Technology
For excellent positioning, the TS offers 45° of tilt technology. The Foot Release Tilt Actuator is easy to operate and prevents tangles caused by exposed cables.

Ability to Grow
With up to 2" of width adjustments and up to 5" of depth adjustments, the TS is able to accommodate the child's growth.



Folding and Rigid Frame Options

The TS is available with a rigid frame for efficiency during propulsion or with a folding frame for easy transport and storage in tight spaces.

Z-finity® System Footplate
Rotating hardware allows infinite adjustment along its rotational arc, achieving multiple angles and placement options.



55° of Rotation-in-Space Technology
For superior positioning, the shortest possible wheelbase, and a smooth tilting motion.



XLOCK® Folding Crossbrace (optional)
Folds compactly for convenient transport and locks when unfolded for superior rigidity.



MONO Backrest System™
The optional MONO Backrest System allows for independent selection of the backrest width and chair width, eliminates interference between backrest hardware and laterals, offers up to 2" of midline backrest offset, and features 35° of angle adjustment. The dynamic function allows for controlled movement of the backrest.

ROTATION-IN-SPACE WHEELCHAIR COMES IN EITHER A RIGID OR FOLDING FRAME AND OFFERS A SUPERIOR COMBINATION OF POSITIONING, MANEUVERABILITY AND PORTABILITY



Foot Release
Eliminate exposed cables by selecting the easy-to-operate Foot Release Tilt Actuator.

ZIPPIE IRIS





ZIPPIE GS



AVAILABLE WITH
FOLDING OR
RIGID FRAMES
AND STANDARD
OR REVERSE WHEEL
CONFIGURATIONS



Reverse Wheel Configuration

To allow children of all sizes to reach the wheels for maximum independence, the GS's main wheels can be placed in standard (rear of the frame) or reverse (front of the frame) configuration. The axle plate may be adjusted vertically and horizontally to allow for additional growth.



Built-in Growth System

The GS easily grows as the child grows. The frame and versatile growing seat pan allow for up to 2" of width growth and up to 4" of depth growth!



Folding or Rigid Frame

Choose the folding frame option for easy transport and storage in tight spaces. Choose the rigid frame option for greater propulsion efficiency.



Special Feature

Optional One-Arm Drive allows maneuverability to be controlled on one side for children with limited function.



Z-finity® System

The Z-finity System's footplate can be adjusted between 0.5" above the seat pan to 21" below the seat pan to provide proper support as children grow. Its rotating hardware allows infinite adjustment along its rotational arc, achieving multiple angles and placement options for positioning the feet, ankles, and knees.



XLOCK™ Technology

Patented locking crossbrace reduces frame play for greater efficiency and less fatigue. XLOCK features one-step push-to-lock and pull-to-fold functions.

CONVENIENCE
OF A FOLDING
WHEELCHAIR
WITH THE
PERFORMANCE
OF A RIGID
FRAME
THANKS TO
XLOCK™
TECHNOLOGY



VersaRail™ Seat Rail

VersaRail offers shared placement of components and accessories, such as seating, armrests, wheel locks, side guards, etc.

Special Feature

2" of width growth and 3" of depth growth built into frame

ZIPPIE X'CAPE



ZIPPIE ZONE



DESIGNED
FOR ACTIVE
CHILDREN,
THIS RIGID
WHEELCHAIR IS OUR
LIGHTEST
ZIPPIE YET!



Our Lightest Pediatric Wheelchair Ever! Using 7000 series aerospace aluminum and Shapeloc technology, the ZIPPIE Zone is incredibly strong and weighs as little as 14 lbs!



Built-in Growth 3" of frame depth adjustment and 2" of seat sling depth adjustment. Frame width growth is as simple as replacing 4 tubes!



Child-Specific Frame Flared or inset frame and inset backrest fit small children and provide better access to the wheels.

Special Features: Flip-up, angle-adjustable footrest for easy transfers

Push-button hub lock located on seat rail stays out of the way during propulsion

Dynamic Center Rear Caster for greater maneuverability and rearward stability



Small Turning Radius Your child can move shoulder to shoulder with family or friends while driving 3 MPH (walking speed). A smooth riding drive suspension is always active, no matter what seat elevation you drive from. Your child will be able to steal cookies from cabinets, write on the chalkboard, and give you smooches at your level.



Manual Tilt and Recline 50° of Manual Tilt and/or 40° of Manual Recline allow the caregiver to change the child's position or posture to assist with feeding, digestion, respiratory function, and visual orientation.

IMAGINE THAT
THE SKY IS
THE LIMIT
WITH THE
PERFORMANCE,
ULTRA-COMPACT
BASE
AND SMALL
TURNING RADIUS



Powered Seat Elevator 12" of Powered Seat Elevation assists with accessing tall shelves or light switches and interacting with peers at eye level.

ZIPPIE MINI





ARM SUPPORT OPTIONS



T-Post

- Standard and tall posts
- Three height settings



Height-Adjustable

- Peds, kids, adult sizes
- Desk and full length



Locking, Flip-Back

- Locking flip-back operation
- Multiple height positions
- Length-adjustable option



T-Post with Sideguard (available on ZIPPIE X'CAPE)

- Desk and full length classic or foam armpads
- Two height options

Other options, not pictured, are also available.

FOOTPLATE OPTIONS



Kids Angle-Adjustable

- Rigid aluminum plate
- Angle-adjustable



Adult Angle-Adjustable

- Oversized slotted aluminum plate
- Angle-adjustable



90° Footboard

- Removable
- Depth-adjustable
- Angle-adjustable



Flip-Up Platform

- Angle-adjustable
- For use with standard swing-aways



Multi-Angle Adjustable (available on ZIPPIE X'CAPE)

- Compatible with Z-finity™ System
- Depth-adjustable
- Angle-adjustable

Other options, not pictured, are also available.

BACKPOST & FOOT SUPPORT HANGER OPTIONS



Stroller Handles
21" 23" 25"
LOW MED HIGH



Swing-Away, Adjustable Stroller Handles

- Adjusts in height
- Folds in for storage



Removable Extension

- Comfortable hand positions
- Rigidizes back canes
- Mounts upward or downward



Height-Adjustable, Angle-Adjustable Push Handle

- 260° of rotation
- User-friendly thumb button release



Swing-In/Out Extension Tubes
(available on ZIPPIE IRIS & ZIPPIE X'CAPE)

- 80°, 70° & 60° angles
- 3 extension tube lengths



Elevating/Articulating Legrest
(available on ZIPPIE IRIS)

- 2" of height adjustment aligns pivot point with knee joint

Other options, not pictured, are also available.



JAY Zip® Back



One-step release



Back shell color options

Versatile One-Step Release Hardware

- Easy one-step release
- Fits on wide range of tubing: 3/4" with inserts, 7/8" with inserts or 1" without inserts
- One-tool adjustments

JAY Zip® Cushion



Dual-Layered Foam



Cushion cover colors available

Comfortable

- Dual-layered, contoured foam base
- Reticulated foam in the outer cover promotes air movement
- Soft, stretchable outer cover

ZIPPIE Frame Colors




IMPORTANT NOTE:
Print processes cannot accurately represent paint finishes. Actual colors may vary.



*Zebra Kolorfusion™ is available on any frame color. Shown on Candy Blue.


ZIPPIE Q300 M Mini Available Shroud Colors

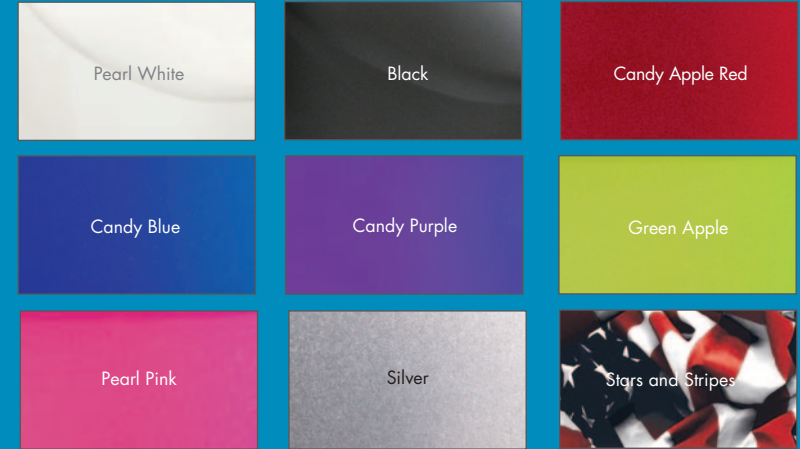

ZIPPIE Q300 M Mini Available Accent Colors



Anodized Parts

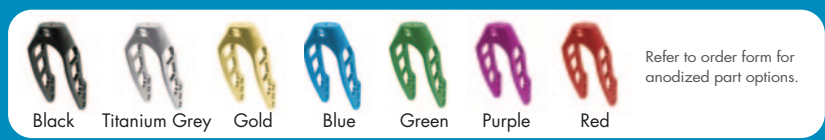

Manual Wheelchair Web Page



Power Wheelchair Web Page



Sunrise Medical (US) LLC, Fresno, CA 93727
932062 Rev.T ©06.2022

800.333.4000 www.SunriseMedical.com

