

INTUITIVE,
EASY-TO-USE
Z | P P | E
VOYAGE
MAKES NO
COMPROMISES
WHEN IT COMES
TO WHAT'S
BEST FOR BABY













### 1. Canopy with Mesh Window

Keep baby shaded during the day, and open the mesh window to check in.

### 2. Angle Adjustable Stroller Handle

One-handed steering that adjusts to caregiver heights.

### 3. Quick Adjust Tilt

Effortlessly tilt forward or rearward for positioning needs.

### 4. Reversible Seating

Easily place the seat facing towards the caregiver, or let baby face outward to interact with their environment.

### 5. Transit Approved

The Voyage has been crash tested and is approved for occupied transit in a motor vehicle.

### 6. Carry-All Storage Basket

Store up to 15 lbs. of essentials.

### 7. Foot Operated Wheel Locks

Hands free! Lock wheels in place easily with your foot.

### 8. Airless Wheels

Flat tires are a thing of the past with airless wheels.

# 9. Swivel Front Wheels with Suspension and Lock-Out

Maneuver around obstacles with confidence, knowing baby is comfortable and secure. Lock-out the front wheels for smooth tracking during longer strolls.





WHEN IT COMES TO SEATING, OUR NUMBER ONE PRIORITY IS POSITIONING

HEALTHY GROWTH

YOUR CHILD FOR

AND DEVELOPMENT



# Voyage Moderate Seating (VMS)

Designed for children with mild to moderate positioning needs. Offers soft, adjustable components and polka dot upholstery with colored accents.



# Voyage Advanced Seating (VAS)

Designed for children with moderate to advanced positioning needs. Offers an extensive selection of seat options and adjustable components, and black stretch upholstery with colored piping.



### Reversible Seat

Face the seat toward you to interact with your child, or forward so your child can explore their surroundings.

### Positioning Accessories

Whitmyer® head supports, lateral trunk supports, anterior trunk supports, medial supports, and hip supports create a seating system that is entirely unique to your child.

Positioning accessories can be moved and adjusted to maintain proper support for your child as they grow.







### Rearward (Posterior) Tilt & Recline

Improves positioning to assist with feeding, digestion, respiratory function, and visual orientation. Posterior tilt range of  $0^{\circ}$  to  $-30^{\circ}$ .

### Forward (Anterior) Tilt

Assists with positioning to encourage healthy development such as bringing hands to the mouth, reaching for objects, or engaging in activities. Anterior tilt range of  $0^{\circ}$  to  $+5^{\circ}$ .

### Height & Angle Adjustable Footrest

Grows with your child and adjusts for various needs throughout the day.

## JCM Hi-Lo Base

Allows your little one to engage in activities at the table or on the floor, while being supported to meet his/her needs.









# Voyage Moderate Seating Colors







Gray







Blue



Red

# FUNCTIONAL ACCESSORIES

The Zippie Voyage is available with options designed to make daily tasks a breeze.



JCM Hi-Lo Base



Folding Booster Base
Height Adjustable



Upper Extremity Support Tray



Utility Hool



Grab Bar



Cup Holde

# Voyage Advanced Seating Colors







Purple



Pink Sorbet



Green



Navy Blue





Turquoise



Black

Butterflies












Battery, Vent or Utility Tray





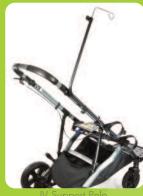

# voyage®



### Technical Specifications

Frame Style:

Early Intervention Device

HCPCS Code: E1232, or E1234 (without seating)

Chair Weight: 28 lbs (base only)

Frame Material: Aluminum User Weight Cap.: 75 lbs Transit Approved: Yes

### Options & Accessories

Frame Type: Tilt, Tilt & Recline

Caster Wheels: 7"

Rear Wheels: 12" Polyurethane

Moderate Seating (VMS): One option

Advanced Seating (VAS): **Backrest Types** 

T-Back Cushion 1" thick T-Back Cushion 2.5" thick

Seating Types Standard Planar with or without hip blocks Contoured Wedge Contoured Anti-Thrust Adjustable Contour

### Frame & Features

(dependent on configuration)

Overall Width: 24"

Seat Widths:

Advanced: 6" to 13" Moderate: 9.5" to 12"

Seat Depths:

Advanced: 5" to 13" Moderate: 8" to 13"

Back Height:

Advanced: (Short) 7" to 13" (Tall) 10" to 17"

Moderate: 16" to 24" (includes head support) Folded Height: 12"

Tilt Range: 35° (5 to -30°)



