


- 
 The Jay[®] Solution[™] cushion that combines the design characteristics from the Jay[®] Easy and Jay[®] Triad cushion lines. It is primarily developed for individuals who require moderate positioning and pressure relief.



MADE IN CANADA



Pre-contouring



Provides pressure reduction, positioning and stability

Positioning Capability



Provides proper positioning for postural needs

Stability



Firmness Layering



Redistributes load while increasing stability and comfort

JAY SOLUTION™

Posture


- The Jay Solution cushion combines the Jay Fluid technology with a molded foam base and visco-elastic foam top layer to optimize the user's comfort and stability.
- The Jay Solution cushion provides dynamic contouring properties and excellent pressure distribution using an inner core of Jay fluid technology.
- The Jay Solution provides enhanced postural pelvic stability with modular pelvic inserts included in the contoured design.
- The Jay Solution has a thigh-contour design to provide the user's legs with gentle adduction support.

Skin

- The Jay fluid pad provides a dynamic contour for the ischials, reducing peak pressures. This comfort and skin protection for users who are at moderate risk for skin breakdown.
- Combine the contoured molded foam base and modular pelvic inserts helps to reduce the potential for sliding and shearing.
- The Jay Solution features a sacral relief designs to reduce pressure at the sacrum for users with a posterior pelvic tilt.

Function

- The use of a stable contoured foam base encourages upper trunk and extremity control.
- Maximizing comfort and positioning, the Jay Solution supports the pelvis to increase user function and increase sitting tolerance.
- The user friendly, moisture resistant inner sleeve ensures cushion durability.
- Maximum weight capacity 250lbs

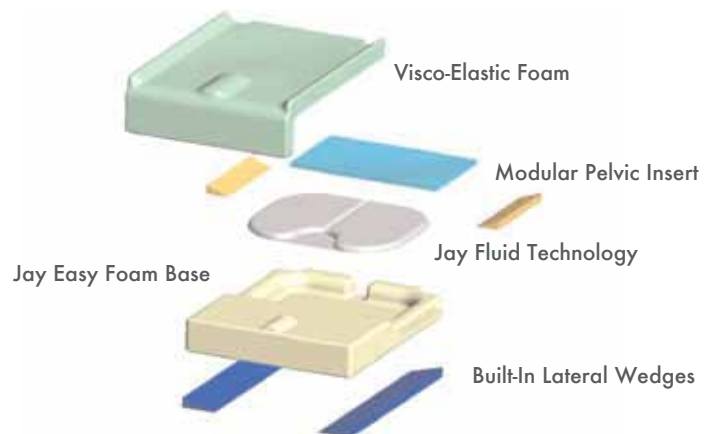

SPECIFICATIONS

JAY®

Funding	ADP: SESJA0360/0365 Under the AADL Program
Weight	4 lbs. (16 x 16)
Product Width	13" - 25"
Product Depth	13" - 22"* The Jay Solution is available in modified format and with solid drop seat.
Product Height	5" 1/4 (at lateral Edge)
Weight Capacity	250 lbs.
Covers	Incontinent

* Certain configurations are required for modification

101733 Rev.B