

Q300M Mini

Compact Size BIG Performance

Q300M Mini

The narrowest
TRUE mid-wheel drive
power wheelchair.

ULTRA-RELIABLE,
HIGH TORQUE 4-POLE MOTORS
WITH SPEEDS UP TO 6 MPH.

OCCUPIED TRANSIT
TESTED TO ISO-7176

SOLID, SINGLE SHEET
METAL CHASSIS

ULTRA LOW 16.5" STFH

SIX WHEEL DRIVE SUSPENSION

CLIMB CURBS
UP TO 3"

RELIABLE

With a chassis formed from a single sheet of metal, steel components with automotive undercoat plating, industry best motors and end of line testing, every Q300 M is built to keep you moving.

TINY TURNING

Thanks to its TRUE mid-wheel drive technology, the Q300 M Mini turns on the spot, for an ultra-small 17" turning radius. Other compact MWD chairs have offset drive wheels, increasing their overall turning radius.

ULTRA-SLIM

Ridiculously narrow at 20.5" wide, the Q300 M Mini goes where others can't. Whether it's a complicated living room, a cluttered kitchen, or a narrow doorway, it'll fit into your life, not the other way around.

ELEVATE AND DRIVE WITH C-ME®

C-Me enables you to elevate up to 12" while traveling up to 3 mph! You can lift or descend at the touch of a button with no need to stop while transitioning. Our suspension doesn't lock out when you drive with C-Me, giving you a smooth and comfortable ride from start to finish.






SEDEO® PRO™

Truly modular in design, the SEDEO PRO seat frame provides a new level of functionality, comfort and design. For moderate to complex clinical needs, the PRO combines a highly adjustable seat frame with an excellent choice of postural support options.



SEDEO® LITE™

The SEDEO LITE is an adjustable seat frame with a sleek look. Its minimal design alleviates weight capacity concerns when using aftermarket equipment. Traditional back canes are width-adjustable from 15" to 20" and easily fit JAY® backs. An array of options cover mild to moderate clinical needs, and intuitive adjustability meets individual requirements.

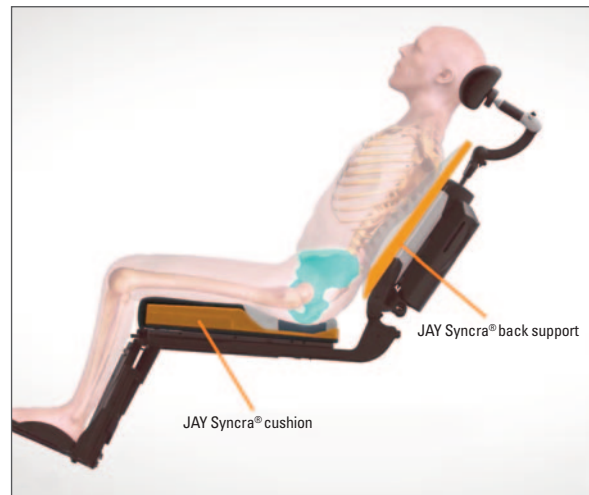

JAY SYNCRA®

Exclusive to the SEDEO PRO Series, JAY Syncra® seat and backrest cushions are specifically designed for optimal performance when tilt is used together with power recline. Featuring a lower rear pelvic wall, deeper leg troughs, a scalloped lower backrest shell and tapered backrest foam, JAY SYNCRA is uniquely engineered to maintain the proper pelvic position throughout the seat movements.



CAPTAIN'S SEAT

The Captain's Seat provides a relaxed sitting position and features a semi-recline, fold-down back, adjustable seat depth, and convertible seat pan that allows you to use your preferred cushion. Available with a 30° tilt and power centermount leg support.

Q300 M SEAT FRAME COMPARISON

	SEDEO PRO	SEDEO LITE	CAPTAIN SEAT
Cryo™ Fluid	Available on JAY Syncra Seating	Available on JAY Fusion	Available on JAY Fusion
Seat Width	16" to 22"	15" to 20"	16", 18", 20"
Seat Depth	14" to 22"	16" to 22"	16" to 22"
Back Support Height	21" to 30"	18" to 30"	25"
Powered Recline	150°		
Anti-shear Back Support	4" (optional)		
Tilt Angle	50°	50°	30°
Elevate with C-Me	12"	12"	
JAY & WHITMYER	●	●	●



Additional information, videos and selection guide on the seating options.

Q300 M Mini			
See order form as specifications may vary based upon configuration			
HPC codes	No Power Seating Options	Single Power Seating Options	Multi Power Seating Options
Group 3 Standard Duty	K0848 / K0849	K0856 / K0857	K0861
Max. user weight	300 lbs		
Seat-to-floor height	16.5" to 19"		
Total width	20.5"*		
Total length	33"		
Min turning radius	17" (Excludes footplates)		
Max. speed	6 mph		
Battery size	Group 22		
Overall weight	Starting at 250 lbs	Starting at 274 lbs	Starting at 322 lbs
Drive wheel position	Mid		
Ground clearance	3"		
Wheel options	12" drive wheels, 14" drive wheels optional; 7" casters		
Max. range	Up to 17 miles		
Transit Brackets	ISO 7176-19:2008 Transit Approved		
Warranty	Lifetime Limited		

* Exact measurements may vary depending on options.

COLOR OPTIONS

SHROUD COLORS



COLOR ACCENTS



Important Note—
Print processes cannot accurately represent painted colors. Actual colors may vary.



Order Form



Web Page



Sunrise Medical (US) LLC, Fresno, CA 93727
MK-129903 Rev.D ©09.2022

800.333.4000 www.SunriseMedical.com

